



BIOCURVE FLOW



EXTRA FLOW SILICONE / short term



Curved
(after insertion)


FEATURES & BENEFITS

- + **Soft Silicone material** – White, radio-opaque, high biocompatible, long term durability and soft tip to avoid vessel wall irritation.
- + **Curved** – If wished can be curved due to silicone soft material. After insertion catheter takes on a curved form. Adding to patients comfort and feel good factor.
- + **Can be inserted over the wire (Seldinger Technique)** – The catheter also includes a removeable stylet from the venous side to allow an easy insertion of the catheter.
- + **Silicone connection tubes** – Durability of catheter.
- + **Highly kink resistant oval shaped Silicone tube** – Excellent constant blood flows.
- + **Excellent flows** – The 15,5 Fr. comparable French size ensures excellent flow rates due to maximum lumen cross section area.
- + **Tip design** – Optimal placing of venous and arterial in/outlets minimizes recirculation – **available with or without side holes.**
- + **No side holes** – Ensures for those using lock solutions an excellent through flow from connector to distal tip. Guaranteeing an excellent filling of catheters between dialysis.
- + **PEEK™ Luer-Lock connector** – High resistance to chemicals and compared to standard including Ultem Luer-Lock connectors even more break resistant.

Catheter tip



D-Line Perfection

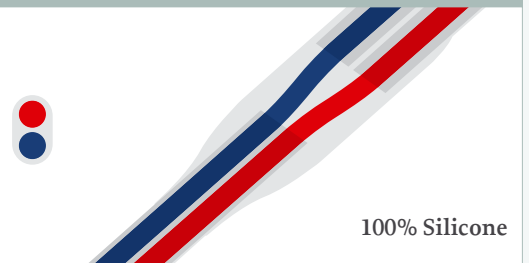


100% Silicone

Surface area within catheter ensures smooth and continuous laminar blood flow, i.e. high flow rates with minimum risk of thrombosis.

FLOW FAMILY ST

NEW

FLOW-COMPARISON

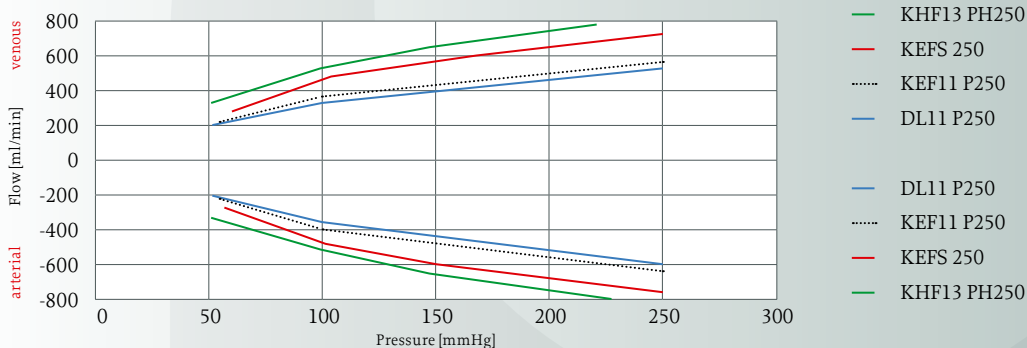

ORDER INFORMATION

BIOCURVE FLOW

EXTRA FLOW Silicone

Item-No.	Outer- ϕ [mm]	Length [mm]	Qty.
N			
KEFS150	• 6,1 x 3,5 OVAL	150	3
KEFS200	• 6,1 x 3,5 OVAL	200	3
KEFS250	• 6,1 x 3,5 OVAL	250	3

N = no sideholes

The Catheter Kit includes:

- 1 x Catheter with stylet
- 1 x Dilator
- 1 x Split Sheath Dilator
- 1 x J-Guide Wire – GWJ 97/700
- 1 x Introducer Needle (IN 18/70)

Quality first:

Joline is DIN EN ISO 9001:2000 and 13485:2003 certified.
 All Joline-Products are CE-Marked.
 Subject to alterations (51852/400002/04/09_RV000)

JOLINE GmbH & Co. KG
 Neue Rottenburger Strasse 50
 D-72379 Hechingen
 Phone +49 (0) 74 71- 98 81- 0
 Fax +49 (0) 74 71- 98 81-111
 www.joline.net, info@joline.net

CE 0124

DISTRIBUTOR